

## At-Home Kit Process for Users: Step By Step

1. Login to Aura in smartphone application or through web app at <https://app.auratracker.org>
2. Once logged in, click on **My Test Kits** tab and click **Add a Test Kit**
  - a. If you already have a kit that looks like 10039 below, click **Activate**



3. Enter **Activation Code** (Barcode Number on Collection Tube)
  - a. 2 letters and 9 digits
    - i. **Label your collection tube with your first name, last name and date of birth**
4. Enter **FedEx Return Tracking Number (12 digits)**
5. Select either **This kit was sent directly to my home** or **This kit was provided to me by my employer**
6. Click **Activate Kit and Begin Collecting**

The screenshot shows the activation form in the Aura app. At the top is a button with a barcode icon and the text 'Scan Bar Code on Tube'. Below this is a horizontal line with the word 'or' in the center. Underneath are four input fields, each with a label and a red asterisk indicating a required field. The first field is labeled 'Activation Code (e.g. Gx12345678)\*'. The second field is labeled 'Barcode Number on Collection Tube'. The third field is labeled 'Collection Date \*' and contains the text '11/19/2021 15:20'. The fourth field is labeled 'FedEx Return Tracking Number (12 digits)\*' and contains the text 'On FedEx return label following TRK#'. Below the input fields are two radio button options. The first option is 'This kit was sent directly to my home'. The second option is 'This kit was provided to me by my employer', which is selected with a blue dot. At the bottom of the form is a large orange button with the text 'Activate Kit and Begin Collecting'.

## FRONT FLAP

# Your GENETWORx COVID-19 Nasal Swab Test Kit Instructions

Read all directions thoroughly and carefully before you start.

- » Before starting Step 1, you must go to [genetworx.auratracker.org/activation](https://genetworx.auratracker.org/activation), or scan the QR code on the right to activate your test kit. The lab cannot process your sample if this step is missed.
- » Once on the site, you will log into your account and activate your kit by entering the barcode number located on the collection tube. After activating your test kit, an instructional video will be available for you to view before starting your collection.



**WARNING:**

1. Do not touch or drink the liquid in the collection tube.
2. In case of spills/splashes into eyes or skin, immediately wash with lots of water. If any symptoms arise after exposure, call a physician.

## SECOND FLAP AFTER OPENING FRONT FLAP

### What's inside the kit:



## BACK SIDE



4060 Innslake Drive, Glen Allen, VA 23060

For questions or concerns regarding your test kit, visit [genetworx.com/contact](https://genetworx.com/contact) or call **610-726-1205**.

For In Vitro Diagnostic Use Only. For Emergency Use Authorization Only. Rx Only. For Use by Individuals 18 Years or Older.








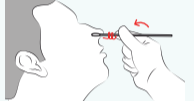








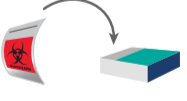

This product (home collection kit in combination with the authorized test) has not been FDA cleared or approved, but has been authorized for emergency use by FDA under an EUA for use by authorized laboratories.

This product has been authorized only for the detection of nucleic acid from SARS-COV-2, not for any other viruses or pathogens.

The emergency use of this product is only authorized for the duration of the declaration that circumstances exist justifying the authorization of emergency use of in vitro diagnostics for detection and/or diagnosis of COVID-19 under Section 564(b)(1) of the Federal, Food, Drug, and Cosmetic Act, 21 U.S.C. § 360bbb-3(b)(1), unless the declaration is terminated or authorization is revoked sooner.

## INSIDE

**Read all instructions and plan same-day Fedex pickup before you start.  
Follow steps 1-18 in order.**

- 1**  Wash and dry your hands before handling the contents of the collection kit.
- 2**  Using a non-smudging blue or black pen/marker, write your first and last name and date of birth on the collection tube.
- 3**  Next, using the same non-smudging pen/marker, fill out the Specimen Confirmation Form.
- 4**  Screw off the top of the collection tube. Place the cap with the opening face up on a clean surface. Do not drink the liquid.
- 5**  Take the swab out of its package. Do not touch the tip of the swab with your hands, and do not lay it down on any surface.
- 6**  Hold swab in one hand and collection tube in the other, being careful not to spill the liquid.
- 7**  Insert the tip of the swab into one nostril, just until the tip of the swab is no longer visible.
- 8**  Rotate the swab while gently circling and pressing around the entire inside edge of your nostril at least 5 times. When finished, take the swab out of your nostril.
- 9**  Using the same end of the swab, repeat step 8 in your other nostril. When finished, take the swab out of your nostril.
- 10**  Place the swab in the collection tube, so that it is submerged in the liquid inside the tube.
- 11**  Break the swab along the pre-scored line by bending it against the collection tube.
- 12**  Screw the top of the collection tube back on, ensuring that it is closed tightly.
- 13**  Discard the broken top part of the swab.
- 14**  Wash and dry your hands thoroughly again.
- 15**  Place the collection tube into the absorbent pad inside the biohazard specimen bag.
- 16**  Seal the biohazard specimen bag by closing the zip lock seal.
- 17**  Fold the biohazard specimen bag so that it fits into the GENETWORx box that this at-home kit came in. Place it inside the box, along with your Specimen Confirmation Form.
- 18**  Place the shipping box into the FedEx return pack, remove the adhesive cover strip and seal the FedEx return pack closed.

### RETURNING THE TEST KIT

**Your test kit must be stored at room temperature and sent back within 24 hours of sample collection.**

Please drop your package off before the last pickup of the day. Do not ship your package late Friday evening, or over the weekend. Do not take it to a FedEx office.

The pre-addressed FedEx return pack can be mailed to GENETWORx via calling FedEx to arrange a pick up or by dropping it in a FedEx drop box.

Visit [bit.ly/fedexdropbox](https://bit.ly/fedexdropbox) online to view FedEx drop box locations and pickup schedules.